



AMAZING SOUND
EASY RECHARGEABILITY



WIDEX MAGNIFY™

WIDEX

SOUND LIKE NO OTHER



A TRUER FITTING FOR ANY EAR-TIP

with **TruAcoustics™**

For years, hearing aid fittings have relied on predictive models to estimate vent leakage. But the introduction of instant ear-tips and more open fittings required a more sophisticated method that took into consideration the behavior of different ear-tips and how they respond in terms of insertion loss and vent effect.

TruAcoustics is an intelligent algorithm by which the in-situ acoustics of the individual user's ear canal can be integrated in the parameter settings of the hearing aid, so that all processes and gain settings are optimized at first fit. **TruAcoustics** calculates fittings more accurately than ever before and delivers perfect Widex processing. **TruAcoustics** enables us to take another leap towards even more natural and amazing sound.



MEET THE MARKET'S SMALLEST RECHARGEABLE RIC SOLUTION

WIDEX **MAGNIFY** is our new product family including our top model mRIC R D which is the smallest rechargeable lithium-ion receiver-in-canal hearing aid on the market.

WIDEX **MAGNIFY** mRIC R D makes recharging simple and easy. A wireless charging solution means no battery door, and thereby fewer gaps for moisture to enter through. Also, key areas of the hearing aid have been treated with an advanced water-resistant nano coating to reduce the risk of moisture damage.

This makes mRIC R D our best protected hearing aid ever.

Convenient, wireless charging and intuitive operation. Handy, portable size.

4-hour charge for full day use.
30-minute charge for 4-hour use.

Moisture resistant nano-coating.

Direct streaming with compatible Apple® and Android™ phones




WIDEX **MAGNIFY**

PERSONALIZATION AND CONNECTIVITY

With WIDEX **MAGNIFY**, users have the option to fully control their hearing aid using the WIDEX **MAGNIFY** app on their Apple or Android smartphone. In addition to many other controls, it's how they access the power of SoundSense Learn.*

WIDEX **MAGNIFY** is also equipped with Real Life Insights so you can give even more personalized counseling. By giving you a sense of the listening environments that are important to your patients, the data makes it possible to offer tailored tips on how to improve the performance of their hearing instrument even further. With their consent, you can use Real Life Insights to download your users' preferences into **Compass™ GPS**.

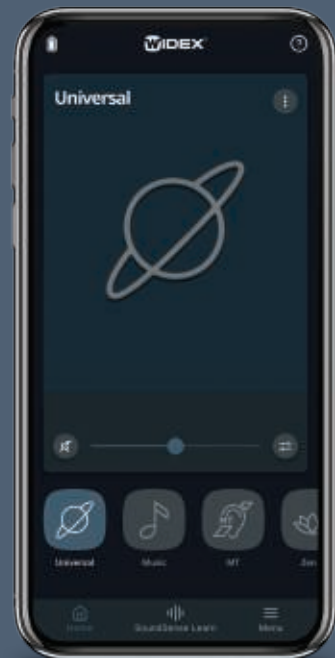

IOS

- ✓ WIDEX **MAGNIFY** app
- ✓ Direct streaming



ANDROID

- ✓ WIDEX **MAGNIFY** app
- ✓ Direct streaming




ARTIFICIAL INTELLIGENCE VERY REAL BENEFITS

WIDEX **MAGNIFY**

With WIDEX **MAGNIFY**, you can now give your users the revolutionary AI-powered **SoundSense Learn**[®] personalization. That means more control over their hearing for your patients, and fewer in-person consultations for you.

To date, **SoundSense Learn**[®] has created almost one million personal sound solutions for users around the world – each one making our AI-engine more efficient than ever. It works so well that 74 percent of users say it improves the listening situation, and 82 percent would recommend it (WIDEXPRESS #41).

*SoundSense Learn is available on the Magnify 440, 330, and 220 models.

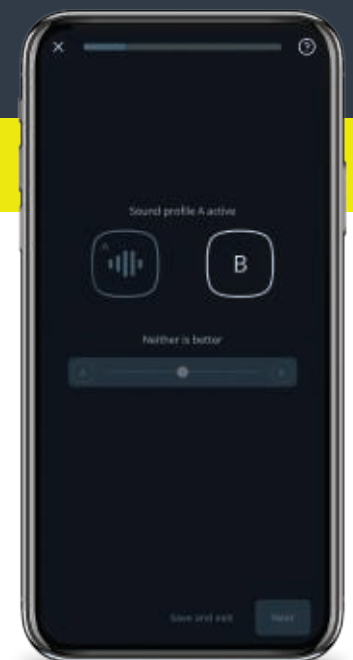



COMPLETE CONNECTIVITY

With WIDEX **MAGNIFY**, connectivity isn't limited to our apps. Our range of communication solutions keeps users connected to a world of sound.



WIDEX TV PLAY™



reddot award winner

Stylish and Stable TV Streaming with WIDEX TV PLAY

WIDEX **MAGNIFY** 2.4 GHz hearing aids work in perfect unison with WIDEX **TV PLAY**, our versatile, stable and straightforward streaming solution.

FEATURES

- Maximum streaming stability lets users move freely around the room without missing a second of the outstanding stereo TV sound that **TV PLAY** delivers
- An elegant and versatile TV sound-streaming solution, **TV PLAY** integrates nicely with modern décor and can be placed in front of the TV or attached to the back
- Intuitive to install and seamless to pair, users can easily customize their listening experience with the app and enjoy their favorite TV moments

You can control **TV PLAY** with the button on your hearing aid, your RC-DEX remote control or the WIDEX **MAGNIFY** app.

COM-DEX

A true hands-free device for streaming high-quality sound from any mobile phone to hearing aids. Can be combined with the COM-DEX app for convenient control of hearing aid volume, sound direction and program settings.

FEATURES

- True hands-free communication
- Compatible with any mobile phone
- Stream audio from all Bluetooth devices
- Available in three stylish colors
- 8 hours of streaming time
- 8 days of standby/control time*
- Comfortable fabric neck-loop
- Control your hearing aids via COM-DEX app

*12 hours per day



COM-DEX REMOTE MIC

COM-DEX Remote Mic is a stylish add-on to your COM-DEX solution and is the answer to difficult one-to-one situations, where clarity of speech is essential.

FEATURES

- Connects wirelessly via COM-DEX to all Widex wireless hearing aids
- Convenient steel clip for easy attaching
- Mute/Unmute microphone directly on the COM-DEX Remote Mic button
- Microphone can also be controlled using the COM-DEX app
- Up to 8 hours of streaming time



RC-DEX

A stylish, compact and user-friendly remote control that gives users better control of selected hearing aid functions. Allows the users to adjust volume and change between programs quickly and discreetly. Small enough to fit on a keyring. A lock switch prevents the remote control being activated accidentally.

FEATURES

- Program toggle
- Volume up/down
- Keyring attachment
- Approximately 12 months of battery life



FM+DEX

A high-quality streaming device that is small and discreet, and designed specifically for wireless hearing aids. Has three functions: FM, telecoil and audio line in, giving the users great flexibility. Easily operated with a toggle button.

FEATURES

- A flexible device for FM streaming
- Small and easy to handle
- Antenna range up to 30 m
- 10 hours of battery capacity



TV-DEX

A palm-sized, user-friendly wireless assistive listening device purpose-built for enjoying TV and audio. Allows up to 10 hours of uninterrupted stereo TV or hi-fi sound between recharges.

FEATURES

- High-quality audio
- Ultra-low delay "Echo-Free™"
- Inputs for both TV and hi-fi systems
- Room Off feature
- 10 hours of non-stop capacity
- Rechargeable using the TV-Base



PHONE-DEX 2

FEATURES

- Integrated answering machine
- Range from base up to 164 feet
- Stores up to 150 contacts
- 320 hours of standby time and 14 hours of call time
- Additional handset available separately



HEARING AID PERFORMANCE & FEATURE OVERVIEW

AUDIBILITY/LISTENING IN QUIET	MAGNIFY 440
TruAcoustics with Personal Gain Integrator	Yes
Digital Pinna	Yes
Enhanced Dynamic Range Compression	Yes
Fitting rationale options	WFR/NAL/DSL
LISTENING IN NOISE	
True-Input Technology	Yes
High Definition Locator (Automatic Adaptive Directionality)	15 ch
Speech Enhancer	Yes (RT/InterEar)
Noise Reduction	Yes
Soft Level Noise Reduction	Yes
Directional Focus	Via MAGNIFY app
Reverse Focus program	Via MAGNIFY app
USER COMFORT	
Preference Control	Yes
Fluid Sound Analyzer (sound classes)	11 (InterEar)
Multi-directional Active Feedback Cancelling	InterEar
TruSound Softener	Yes
SMARTWIND MANAGER™/Wind Noise Attenuation	Yes
USER CONTROLS	
Listening programs	5
Speech (S)/Music (M) streaming template	Yes
Speech (S)/Music (M) streaming template	Yes
Phone+ (available in select styles)	Yes
SmartSpeak/SmartTone	Yes
Zen/Zen+	InterEar
HEARING HIGH FREQUENCIES	
High-Frequency Boost	Yes
Audibility Extender	Yes
FITTING AND FINE TUNING	
Processing channels	15
Fitting frequency channels	15
Sensogram	Yes
Sound Diary/Datalogging	Yes
APPS, ACCESSORIES AND REMOTE CONTROLS	
MAGNIFY APP WITH SOUNDSSENSE LEARN ¹	Yes ¹
TONELINK APP	Yes
TV PLAY ¹	Yes ¹
COM-DEX	Yes
COM-DEX REMOTE-MIC	Yes
RC-DEX	Yes
TV-DEX	Yes
FM+DEX	Yes


¹Only mRIC R d, RIC 312 D, BTE 13 D and BTE R D models ²MAGNIFY APP without SOUNDSSENSE LEARN.

MAGNIFY 330	MAGNIFY 220	MAGNIFY 110
Yes	Yes	Yes
Yes		
Yes	Yes	Yes
WFR/NAL/DSL	WFR/NAL/DSL	WFR/NAL/DSL
Yes	Yes	Yes
12 ch	10 ch	
Yes		
Yes	Yes	Yes
Yes	Yes	Yes
Via MAGNIFY app	Via MAGNIFY app	Via MAGNIFY app
Via MAGNIFY app	Via MAGNIFY app	
Yes	Yes	Yes
7 (InterEar)	4	3
InterEar	InterEar	Yes
Yes	Yes	
4	3	3
Yes	Yes	Yes
Yes	Yes	Yes
Yes	Yes	Yes
InterEar	InterEar	InterEar
Yes	Yes	Yes
12	10	6
12	10	6
Yes	Yes	Yes
Yes	Yes	Yes
Yes ¹	Yes ¹	Yes ²
Yes	Yes	Yes
Yes ¹	Yes ¹	Yes ¹
Yes	Yes	Yes
Yes	Yes	Yes
Yes	Yes	Yes
Yes	Yes	Yes

THE WIDEX MAGNIFY FAMILY




Model	RIC 10		RIC 312 D			
Style type	RIC		RIC/RITE			
Receiver	S-receiver	M-receiver	S-receiver	M-receiver	P-receiver	HP-receiver
Open fit option (BTE/RIC)	Yes	Yes	Yes	Yes		
Tube type						
Battery size	10	10	312	312	312	312
Battery life*	100 hrs	100 hrs	165 hrs	165 hrs	160 hrs	155 hrs
Dust and weather protection	IP-68	IP-68	IP-68	IP-68	IP-68	IP-68
Rechargeable solution						
USER CONTROLS						
Preference control	via DEX or TONELINK App	via DEX or TONELINK App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App
Program shift	via DEX or TONELINK App	via DEX or TONELINK App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App
Telecoil	via FM+DEX	via FM+DEX	Yes	Yes	Yes	Yes
ACCESSORIES						
DEX device compatible ¹	All	All	All	All	All	All
Remote control (optional)	via DEX or TONELINK App	via DEX or TONELINK App	via DEX or MAGNIFY App	via DEX or MAGNIFY App	via DEX or MAGNIFY App	via DEX or MAGNIFY App
TV PLAY			Yes	Yes	Yes	Yes
DAI compatible	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX
TECHNICAL DETAILS						
Max Gain/MPO (dB/dB SPL) ²	63-118	69-124	62-118	69-124	70-132	84-138
Bandwidth ³ (Hz)	100-10000	100-10000	100-10000	100-10000	100-6750	100-7300
Fitting range (dB)	0-80	0-90	0-80	0-90	40-110	50-115

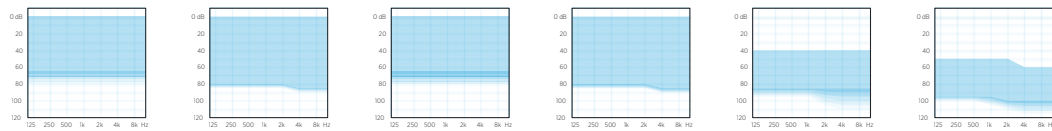

*Battery Life in real-life situations depend on hearing aid features used and streaming time.
¹COM-DEX, COM-DEX REMOTE MIC, RC-DEX, TV-DEX, FM+DEX. ²IEC 60118-0 Peak values. ³According to 60118.0:1983



Model	mRIC R D	BTE R D	BTE 13 D	XP	CUSTOM IM/IP	CIC-M	CIC
Style type	RIC/RITE		BTE	BTE	ITE/ITC	CIC MICRO	CIC
Receiver	V.2 M-receiver	V.2 P-receiver			IM/IP receiver		
Open fit option (BTE/RIC)	Yes	Yes	Yes	Yes			
Tube type			V.2 Thin/Classic	V.2 Thin/Classic			
Battery size	Li ion/312	Li ion/312	Li ion/13	13	312	312	10
Battery life*	40 hrs	35 hrs	50 hrs	235 hrs	150 hrs	IM: 160 hrs IP: 135 hrs	110 hrs
Dust and weather protection	IP-68	IP-68	IP-68	IP-68	IP-68	IP-68	IP-68
Rechargeable solution	Li ion	Li ion	Li ion				
USER CONTROLS							
Preference control	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via DEX or TONELINK App	via button, DEX or TONELINK App	via TONELINK App
Program shift	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via button, DEX or MAGNIFY App	via DEX or TONELINK App	via button, DEX or TONELINK App	via TONELINK App
Telecoil	via FM+DEX	via FM+DEX	Yes	Yes	Yes	via FM+DEX	via FM+DEX
ACCESSORIES							
DEX device compatible ¹	All	All	All	All	All	All	All
Remote control (optional)	via DEX or MAGNIFY App	via DEX or MAGNIFY App	via DEX or MAGNIFY App	via DEX or MAGNIFY App	via DEX or TONELINK App	via DEX or TONELINK App	via TONELINK App
TV PLAY	Yes	Yes	Yes	Yes			
DAI compatible	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX	via FM+DEX
TECHNICAL DETAILS							
Max Gain/MPO (dB/dB SPL) ²	69-128	71-133	73-140	74-140	67-128	IM: 63-123 IP: 68-127	63-124
Bandwidth ³ (Hz)	100-10000	100-8050	100-7900	100-7750	100-9800	100-10000	100-10000
Fitting range (dB)	0-90	20-110	10-105	10-105	0-85	Up to 90	up to 90

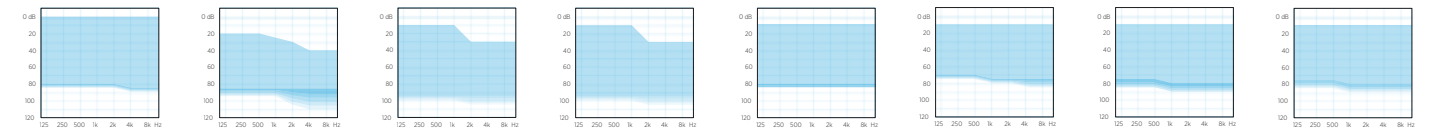

HEARING AID COLORS



MAGNIFY	123	119	121	118	116	115	120	114	068	077	076	136	078	025	107	028
RIC/BTE Models	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•
ITE Models															•	•



WIDEX MAGNIFY



SOUND LIKE NO OTHER



Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play and Google Play logo are trademarks of Google Inc.

MYWIDEXPRO.COM

Subject to change with no further notice.
Widex is a registered trademark of WSAUD A/S.
DBR01MM/0722